AQUABRASION DEVICE
Model Number: BD0030EBCOE

User Manual
Ръководство за потребителя
Užívateľská príručka
Brugsanvisning
Benutzerhandbuch
Gebruiksaanwijzing
Kasutusjuhend
Εγχειρίδιο χρήστη
Manual del usuario
Guide de l'utilisateur
Korisnički priručnik
Manuale Utente
Lietotāja rokasgrāmata
Naudotojo vadovas
Használati útmutató
Gebruiksaanwijzing
Brukerveiledning
Podręcznik użytkownika
Manual do utilizador
Manual de utilizare
Navodila za uporabo
Užívateľská príručka
Käyttöohje
Bruksanvisning
Kullanım Kılavuzu
Посібник користувача

English	3
Български език	14
Čeština	25
Dansk	36
Deutsch	47
Eesti keel	58
Ελληνικά	69
Español.....	80
Français	91
Hrvatski	102
Italiano	113
Latviešu valodā	124
Lietuvių kalba	135
Magyar	146
Nederlands	157
Norsk	168
Polski	179
Português	190
Română	201
Slovenščina	212
Slovenský jazyk	223
Suomi	234
Sverige.....	245
Türkçe.....	256
Українською	267

USER MANUAL

ARTISTRY LABS™ Aquabrasion is intended to be used as part of ARTISTRY LABS Retexturizing System with ARTISTRY LABS Retexturizing Peel and ARTISTRY LABS Retexturizing Serum. See page 10 for complete Treatment System instructions. This 28 Day Skincare System is recommended to be used every quarter.

NOTICE: Before using the device for the first time, it is essential to read the instructions and safety information in this booklet to ensure proper use and to achieve your best results. We recommend saving these instructions for future reference.

IMPORTANT SAFEGUARDS

The following basic safety precautions apply when using electrical appliances, please read all instructions before use. And keep this instruction manual for future use.

- (1) This product is intended for household and indoor use only. Do not use this product for commercial use, and do not use it outdoors.
- (2) This appliance must be used exclusively by adults. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be done by children without supervision.

WARNING

Do not stare at the operating light source.


Table of Contents

- Device Information
- Charging Information
- SUCTION Mode
- SERUM Mode
- Cleaning Method
- Weekly System Regimen
- Safety & Caution
- Technical Information
- Certifications and Marks of Conformity
- Distributor & Manufacturer Contact Information
- Customer Service Information
- Warranty


Device Information

Components

- Device x 1
- Device Cradle x 1
- Water Tank Holder x 1
- SUCTION Nozzle Head x 2 (Large/Small Nozzle)
- SERUM Dome Attachment x 1
- USB-C Cable x 1
- User Manual x 1
- Cleaning Brush x 1

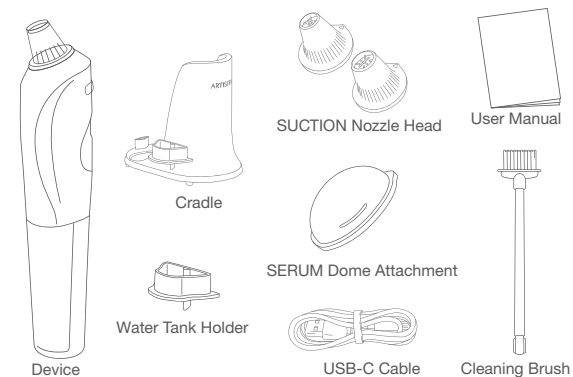Device Structure

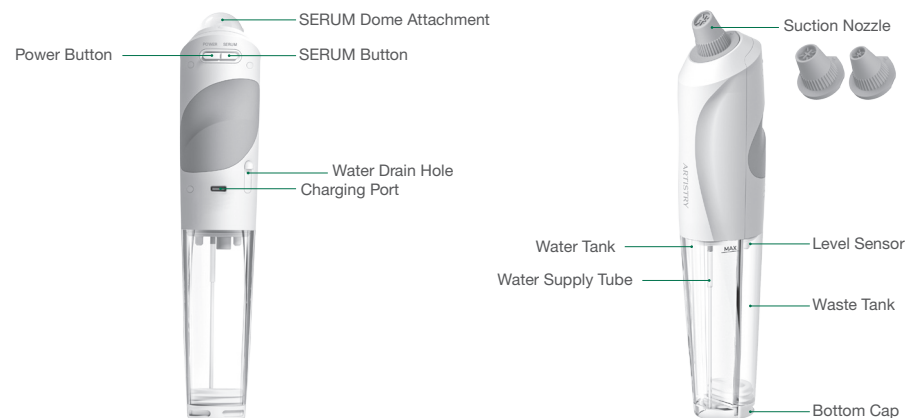Charging Information

Before Use: Charge the Device with USB-C Cable provided. Charging time can vary by country and outlet.

How to Charge

1. Ensure Tanks are empty. Securely place the Device upright in the Cradle.
2. Use the USB-C Cable provided to charge the Device*
3. A red light will be displayed around the POWER Button while charging. When the Device is fully charged, the light will turn green.
4. When charging is complete, unplug it. The Device will not operate while plugged in.

* 5V adapter is required (not included). Do not use a quick/fast USB adapter.

Low Battery Indicator

When the battery Power is less than 5%, the Device will beep several times and a red light will blink around the POWER Button. The Device will not work. Follow the steps above to charge the Device.

Turning the Device On/Off and Changing Modes/Levels

- To turn the Device on, press the POWER Button. A white light will indicate the Device is on.
- To select the appropriate Suction level, press the POWER Button again for Level 1; press twice for Level 2; press three times for Level 3.
- To select SERUM Mode, press the SERUM Button once the Device is on.
- To turn the Device off, press and hold the POWER Button for 2 seconds.

NOTE: You may turn the Device off whenever necessary or desired.

SUCTION Mode Developed to use after ARTISTRY LABS Retexturizing Peel

Suction Nozzle Selection

There are two sizes of Suction Nozzles for different areas of your face.

- Small Nozzle: For smaller crevices like the nose area.
- Large Nozzle: For larger areas like the forehead, cheeks, and chin.

Changing the Suction Nozzle

Ensure the Device is off prior to changing the Suction Nozzles. To remove the Suction Nozzle from the Device, gently pull outward. To replace the Nozzle, align the point on Nozzle with the triangle on Device and firmly push it into place while sliding it onto the ports.

NOTE: When changing Nozzles, it is easier to do when hands are dry.

Suction Levels

The Device has three Suction Levels to choose from depending on your skin tolerance. We suggest starting with Level 1 and then increasing intensity depending on your comfort Suction Level. Turn the Device on and select the appropriate Suction Level by using the POWER Button: Press once more for Level 1; press twice for Level 2; and press three times for Level 3.

LEVEL	DESCRIPTION	INDICATOR LIGHT
Level 1	Low Suction	Blue 
Level 2	Medium Suction	Yellow 
Level 3	High Suction	Pink 

Filling the Water Tank

1. Remove the Bottom Cap.
2. Firmly grip the top half of the Device and pull the Water Tank down with the other hand to remove it from the Device.
3. Fill the Water Tank with water until it reaches the MAX limit line*.

*Do not overfill past the MAX limit line.

4. Push the Water Tank securely back into place.
5. Replace the Bottom Cap.


Suction Treatment Instructions

1. Place the appropriate Suction Nozzle on the Device.
2. Fill the Water Tank with water.
3. Turn the Device on by pressing the POWER Button.
4. Choose the appropriate Suction Level by pressing the POWER Button again accordingly.

*The Device will take a few seconds to “ramp up” and get up to speed.

5. Gently sweep the Device across the face ensuring the Suction Nozzle is in good contact with the skin. Follow the diagram below for directions for treatment.
6. Continue until the entire face has been treated.

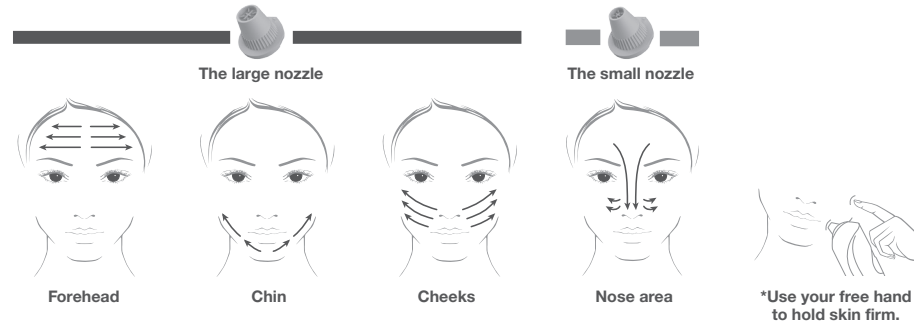NOTE:

- Avoid the eye area.
- Do not hover the Device in a single spot for more than 3 seconds.
- If the Device is tipped too much, the Device will shut off. Restart the Device and continue treatment.
- Stop using it if you feel excessive discomfort.
- You may experience slight redness of your skin during use. This is normal.
- Do not operate in Suction Mode without water in the Water Tank.
- Empty the Waste Tank and refill the clean Water Tank as many times as necessary throughout treatment. The Waste Tank is not removable from the Device.
- Do not shake, tip, or lay the Device down when the Waste Tank is filled with water.
- The Device will automatically turn off after 5 minutes of use in Suction Mode to save power. If needed, turn the Device back on and continue the treatment.


When Suction Treatment is Complete

1. Press the POWER Button for 2 seconds to turn the Device off.
2. Hold the device upright over the sink then slightly tilt to avoid spilling over your hands.
3. Remove the bottom cap from the Waste Tank to allow the wastewater to drain then replace the cap and follow cleaning instructions.
4. Ensure the Device and Suction Nozzles are properly cleaned and dry.

SERUM Mode Developed to use after ARTISTRY LABS Retexturizing Serum

SERUM Treatment Instructions

1. Attach the SERUM Dome Attachment to the Device by aligning the point on the Dome Attachment with the triangle on the Device.
2. Push the Dome Attachment into place while sliding onto the ports.
3. Turn the Device on by pressing the POWER Button.
4. Press the SERUM Button. A red light will appear inside the Dome Attachment.
5. Starting at the top of the forehead, lightly touch the Dome on the face. The Device will beep after 3 seconds. Move 2-3cm (1 inch) and repeat until the entire face is treated.

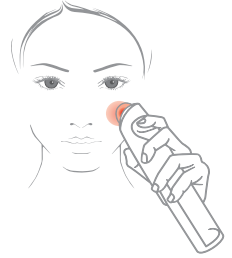NOTE:

- Avoid the eye area.
- Do not shine LED light into your eyes. If necessary, close your eyes during treatment.
- SERUM Mode will last 10 minutes, but you may turn the Device off at any time.

When the SERUM Treatment is Complete

1. Press the POWER Button for 2 seconds to turn off the Device.
2. Grasp the Dome Attachment by the sides with two fingers and gently pull it outward and rinse it to clean.
3. Return the Device to the Cradle and plug it in to charge.

Cleaning Method


Please turn off the Device before cleaning.

This appliance should be cleaned after every use to avoid the accumulation of grease and other residues.

Cleaning and user maintenance shall not be done by children without adult supervision.

Cleaning Mode after Using the SUCTION Mode

- Place the Large Suction Nozzle on the Device.
- Fill the Water Tank with water until it reaches the MAX limit line.
- Turn the Device on by pressing the POWER Button then press and hold the SERUM Button for 2 seconds to start the Cleaning Mode.
- Turn the Dome Attachment upside down and securely hold it on top of the Large Suction Nozzle.
* Hold the Device upright. If the Device is tilted too much, it will beep and shut off. If this occurs, hold the Device upright and turn it on to continue cleaning.
- When the water in the Water Tank is exhausted, press and hold the POWER Button for 2 seconds to turn the Device off.
- Hold the device upright and remove the Bottom Cap to empty the Waste Tank.
- Separate the Water Tank and place it on the Cradle's Water Tank Holder.
- Remove the Large Suction Nozzle and gently rinse the Suction Nozzles, Dome Attachment, and Bottom Cap with water.
- Return the Device to the Cradle and allow all the components to completely dry.

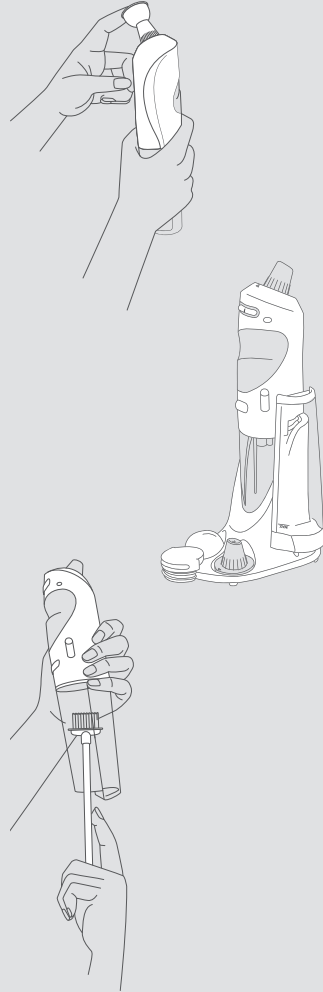How to clean the Waste Tank

- Remove Bottom Cap and empty the Waste Tank.
- Use the brush and dampen it to clean the Waste Tank walls thoroughly and gently.
- Return the Device to the Cradle and allow the Waste Tank to completely dry.

NOTE:

- Avoid soap and excessive/aggressive wiping or brushing. It may cause damage.
- Do not shake, tip, or lay the Device down when the Water Tank or Waste Tank is filled with water. The Device should remain in the upright position if either Tank still contains water.
- Water may seep out of the drain hole. This to protect the pump from water backflow if the Device is tilted too much while the Waste Tank contains water.

[Weekly System Regimen]

ARTISTRY LABS Retexturizing System

This Treatment System consists of two different regimens: Peel Day and Non-Peel Day.

Follow this weekly regimen below for 4 weeks, ensure Peel is used on non-consecutive days:

DAY 1 Peel Day	DAY 2 Non-Peel Day	DAY 3 Non-Peel Day	DAY 4 Peel Day
PEEL DAY INSTRUCTIONS:			
AM:			
<ul style="list-style-type: none"> Apply 2 pumps of Serum to the entire face using fingertips. Let fully absorb and completely dry. 			
PM:			
<ul style="list-style-type: none"> Apply 6-8 pumps of Peel with fingertips onto cleansed skin. Spread as a layer and leave on for 8 minutes. Thoroughly rinse face with water to remove Peel and pat dry. Use the Device + Suction Nozzle to remove skin debris and to extract pores. Pat skin dry. Apply 2 pumps of Serum to the entire face using fingertips. Let fully absorb and completely dry. Follow with Device in SERUM mode + Dome Attachment. 			
DAY 5 Non-Peel Day	DAY 6 Non-Peel Day	DAY 7 Non-Peel Day	
NON-PEEL DAY INSTRUCTIONS:			
AM:			
<ul style="list-style-type: none"> Apply 2 pumps of Serum to the entire face using fingertips. Let fully absorb and completely dry. 			
PM:			
<ul style="list-style-type: none"> Apply 2 pumps of Serum to the entire face using fingertips. Let fully absorb and completely dry. Follow with Device in SERUM mode + Dome Attachment. 			

Skincare products and Device should avoid contact with eyes. If at any time persistent redness, burning, or stinging occurs, discontinue use and/or rinse face. Follow AM and PM steps with a moisturizer.

Order of Use with your daily Skincare Regimen:

- Daytime Regimen: Cleanser, Toner, Retexturizing Serum, Moisturizer and/or SPF.
- Nighttime Regimen with Peel: Cleanser, Retexturizing Peel, Device + Suction Nozzle, Toner, Retexturizing Serum, Device in SERUM Mode + Dome attachment, Moisturizer.
- Nighttime Regimen without Peel: Cleanser, Toner, Retexturizing Serum, Device + SERUM Mode + Dome Attachment, Moisturizer.

Safety & Caution

The manufacturer and distributors are not responsible for Device failure or damages as a result of not following the use instructions and warnings set in this manual.

Safety Statement

1. Do not use the Device if you:
 - Have a facial implant (e.g. metal, silicone)
 - Have eczema, dermatitis, or other skin ailments.
 - Are pregnant or breastfeeding.
2. Consult a doctor:
 - If you had recent facial surgery
 - If you are under medical treatment.
3. Do not use the Device on the following areas:
 - Injured, inflamed, or irritated skin
 - Where surgery was performed
 - Eyes and areas around the eyes
 - Where hair removal treatment has recently been performed (less than 24 hours)
 - Areas with pain perception or sensory disabilities

Caution Statements

1. To avoid potential fire hazards, electric shock or injury, or damage to the Device:
 - Do not expose the Device to extreme temperatures (higher than 45°C (113°F) or lower than 0°C (32°F)), heat sources, direct sunlight, or flammable materials.
 - Do not clean with flammable substances such as alcohol, benzene, or thinner.
 - Do not use the Device if it emits smoke, odor, or unusual noise.
 - Do not place the Device in a microwave or high-pressure container.
 - Do not modify or repair the Device.
 - Do not drop, throw, or apply heavy impact to the Device.
 - Do not submerge the Device in water and avoid contact with a large amount of water.
 - Do not charge the Device in a damp place.
 - Use dry hands with charging cable/electrical outlet.
 - Keep the Device out of the reach of children (under 12 years of age).
2. For optimal Device performance:
 - Power off when not in use.
 - Charge the Device by connecting the included USB-C cable to the 5V adapter to charge. Do not use a quick/fast USB Adapter.
 - Turn the Device off before charging.
 - Fully charge the battery at least every 4 weeks when not in use. This will help maintain the battery lifespan.
 - Turn the Device off prior to changing the suction nozzles/dome. A malfunction may occur if the Device is left on while changing the attachments.
 - Do not turn the Device upside down.
 - Do not use detergent or boiling water to cleanse the Suctions Nozzles / Dome Attachment as it may cause discoloration or corrosion.
 - Do not expose the Device or Components to direct sunlight for an extended period of time as it may cause discoloration or affect product performance.

TECHNICAL INFORMATION

Product Name: ARTISTRY LABS™ AQUABRASION

Model Number: BD0030EBCOE

Rated Voltage/Input: 3.7 V DC/ 1.4 A

Contains: Li-polymer 900 mA Rechargeable Battery

USB-C cable included.

Local charger – not included.

Charging indicator – present

Device size parameters (mm)

Body: 231.5(H)*52(W)*45(L)

Cradle: 108(H)*116.5(W)*76(L)

Device weight: 255g (main body), 70g(cradle)

Materials used:

Device Body Plastic: P.C+ABS white

Waster and Water Tank material: ECOZEN (clean)

Nozzle Tips: Silicone

Cap lens: P.C

Operating temperature: 5°C~25°C (41F - 77F)

Storage conditions: away from direct sunlight, in a dry place, out of children's reach

Storage temperature: at -20°C~ +25°C (-4F - 77F)

Battery

Model name: XK603048

Battery Manufacturer: DONGGUAN LEAN POWER NEW ENERGY TECHNOLOGY CO., LTD
NO.198, Jin Chuan Road 2, Xie Gang Town, Dong guan, Guangdong.

Manufacture Date: Marked separately on battery.

A/S center: 031-701-0951

MADE IN KOREA

DEVICE CERTIFICATIONS / MARKS OF CONFORMITY:

EC DECLARATION OF CONFORMITY

This electronic equipment is in conformity with relevant European harmonization legislation on consumer and environmental safety.

WASTE AND ENVIRONMENT

This electronic must be disposed separately from household waste at its life end. Please dispose of it at your local recycling system. Correct Disposal of this product (Waste Electrical & Electronic Equipment) (Applicable in the European Union and other countries with separate collection systems) The European Directive 2012/19/EU on Waste Facilities Electrical and Electronic (WEEE) This marking shown on the product or its literature, indicates that it should not be disposed with other household wastes at the end of its working life. To prevent possible harm to the environment or human health from uncontrolled waste disposal, please separate this from other types of wastes and recycle it responsibly to promote the sustainable reuse of material resources. Household users should contact either the retailer where they purchased this product or their local authority office, for details of where and how they can take this item for environmentally safe recycling.

WARRANTY:

Please visit Amway website product page in your country for warranty details.

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